

SUPPORT MATTERS!

NEW & IMPROVED!

Therapeutic 5 Zone

LONG TERM CARE SUPPORT FOAM MATTRESS



Committed
TO QUALITY

The perfect introduction to long term care therapeutic foam mattresses, the Therapeutic 5 Zone offers high quality latex-free foam for pressure redistribution. Upgrade from an innerspring for better compliance and durability while improving patient satisfaction.

FEATURES

- Deluxe horizontal, cross-cut foam provides support and comfort over 5 pressure zones
- Removable fluid resistant nylon cover is low shear and vapor permeable
- Non-skid bottom keeps mattress safely in place
- Meets flammability standard CFR 1632 and 1633
- Compression technology creates a 25"W x 9.8"D x 36"H box footprint creating convenience for transport and delivery

WHAT'S NEW AND IMPROVED?



Weight Capacity
3rd Party Verified



Made In
North America



Enviromentally-Friendly
Product Design

Specifications


Model #	15019
Size	36" x 80" x 6"
Weight Capacity	350 lbs.
Warranty	3 years limited foam, 1 year cover

SUPPORT MATTERS!

NEW & IMPROVED!

Gravity 7

LONG TERM CARE PRESSURE REDISTRIBUTION FOAM MATTRESS

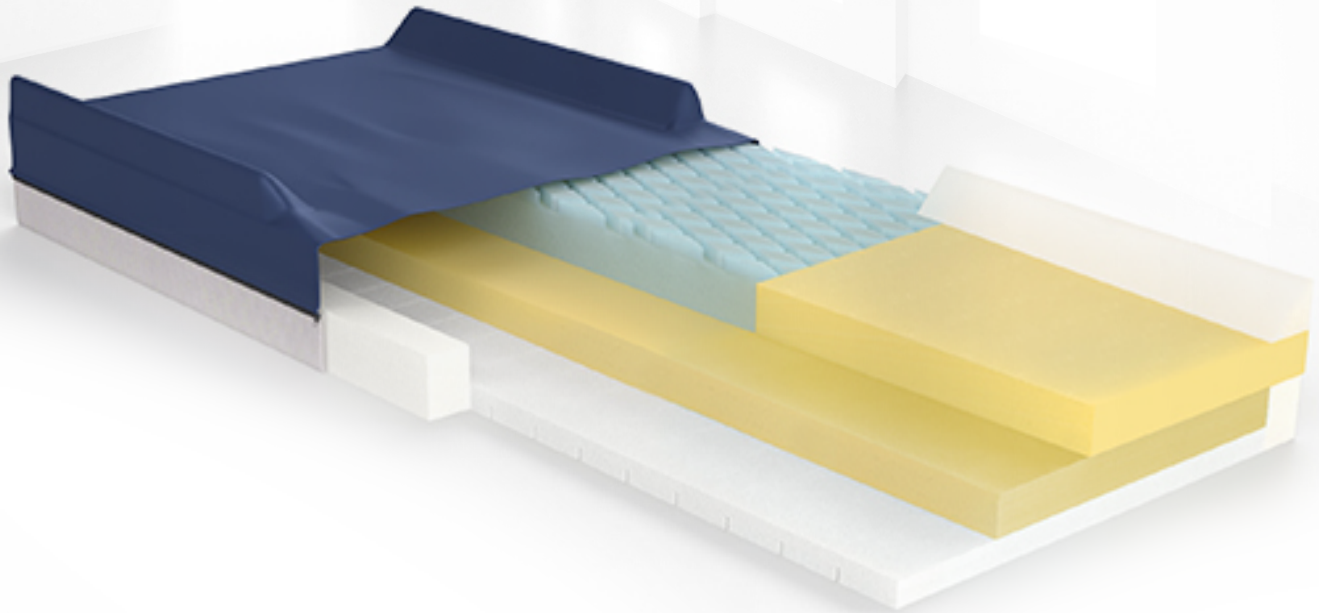

Committed
TO QUALITY

Specifically engineered for pressure redistribution, shear reduction and skin protection, the Gravity series variable density sleep systems include unique combinations of latex-free foam.

FEATURES

- Die-cut full width top layer for pressure redistribution
- Open cell, visco-memory foam heel section for enhanced lower extremity pressure offloading
- Removable premium, multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Available with or without raised rails
- 7 year limited foam warranty, 3 year limited cover warranty
- 15770: Compression technology creates a 25"W x 9.8"D x 36"H box footprint creating convenience for transport and delivery
- 350 lb weight capacity

WHAT'S NEW AND IMPROVED?




3 Year Extended
Cover Warranty*



Weight Capacity
Verified



Made In
North America

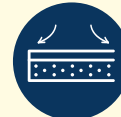

Full Width
Top Layer



15770 Enviromentally-
Friendly Product Design

Model Number	Size	Model Number	Size
15876	36" (W) x 76" (L) x 6" (H)	15777 (with raised rails)	36" (W) x 80" (L) x 6" (H)
15877 (with raised rails)	36" (W) x 76" (L) x 6" (H)	15784	36" (W) x 84" (L) x 6" (H)
15770	36" (W) x 80" (L) x 6" (H)	15785 (with raised rails)	36" (W) x 84" (L) x 6" (H)


* 20,000 cycle static fatigue testing verified

SUPPORT MATTERS!

NEW & IMPROVED!

Gravity 8

LONG TERM CARE PRESSURE REDISTRIBUTION FOAM MATTRESS

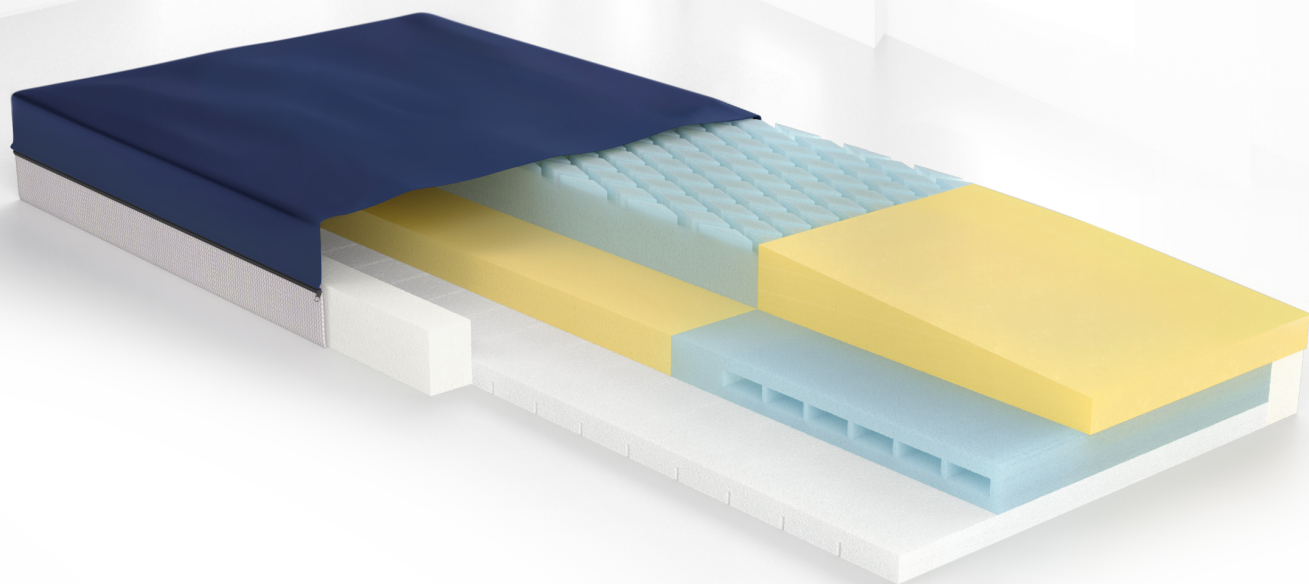

Committed
TO QUALITY

Specifically engineered for pressure redistribution, shear reduction and skin protection, the Gravity series variable density sleep systems include unique combinations of latex-free foam.

FEATURES

- Die-cut full width top layer for pressure redistribution
- Designated memory foam sloped heel section for enhanced lower extremity pressure offloading
- Removable premium, multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Available with raised rails in 80"
- 10 year limited foam warranty, 3 year limited cover warranty
- 400 lb weight capacity

WHAT'S NEW AND IMPROVED?




3 Year Extended
Cover Warranty*



Weight Capacity
Verified



Made In
North America

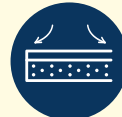

Full Width
Top Layer

Model Number	Size	Model Number	Size
15886	36" (W) X 76" (L) X 6" (H)	15888 (with raised rails)	36" (W) X 80" (L) X 6" (H)
15870	36" (W) X 80" (L) X 6" (H)	15884	36" (W) X 84" (L) X 6" (H)


* 20,000 cycle static fatigue testing verified

SUPPORT MATTERS!

NEW & IMPROVED!

Multi-Ply 6500 Dynamic Elite

LONG TERM CARE PRESSURE REDISTRIBUTION FOAM MATTRESS

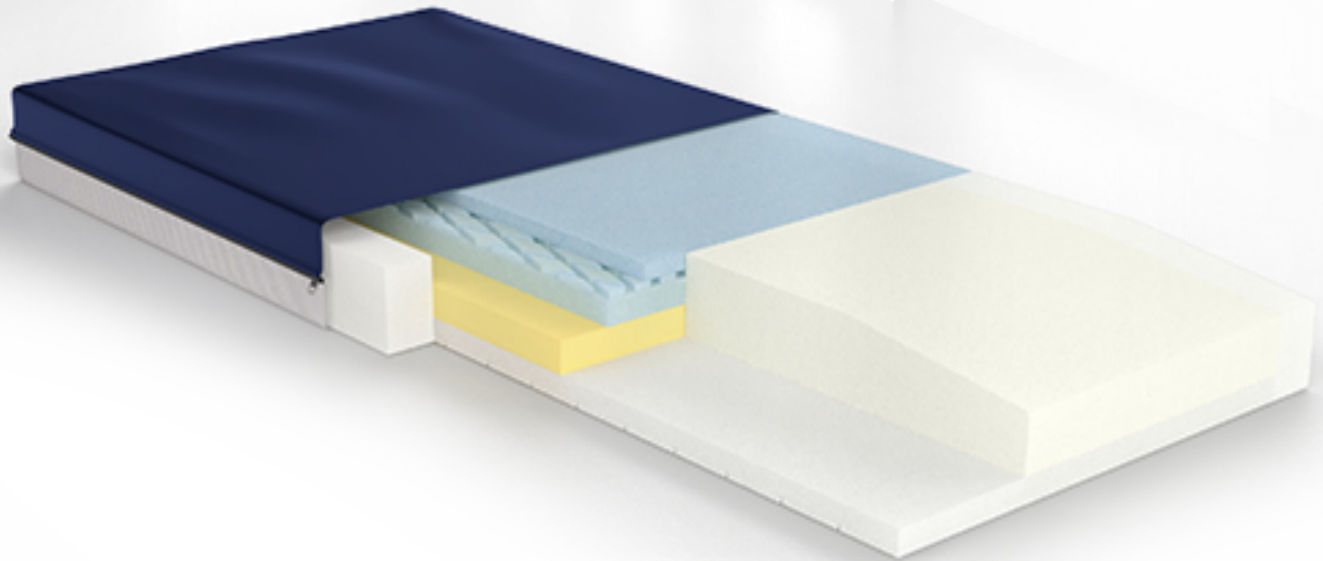

Committed
TO QUALITY

Specifically engineered for pressure redistribution, shear reduction and skin protection, the Dynamic Elite variable density sleep system includes unique combinations of latex-free foam.

FEATURES

- Full width gel infused, open cell visco-elastic foam top layer dissipates heat providing a cool and dry environment to help preserve skin integrity
- Heel slope offloads pressure from vulnerable heels
- Removable premium, multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Full length horizontally scored articulation cut bottom layer extends the durability and life of the mattress
- Meets flammability standard CFR 1632 and 1633

WHAT'S NEW AND IMPROVED?




3 Year Extended
Cover Warranty*



Weight Capacity
Verified



Made In
North America

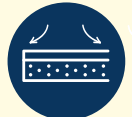

Full Width
Top Layer

Specifications

Model #	6500-DE-2-FB
Size	36" W x 80" (L) x 6" (H)
Weight Capacity	400 lbs.
Warranty	15 years limited foam, 3 year cover


* 20,000 cycle static fatigue testing verified

SUPPORT MATTERS!

NEW & IMPROVED!

Multi-Ply ShearCare 1500

LONG TERM CARE PRESSURE REDISTRIBUTION FOAM MATTRESS

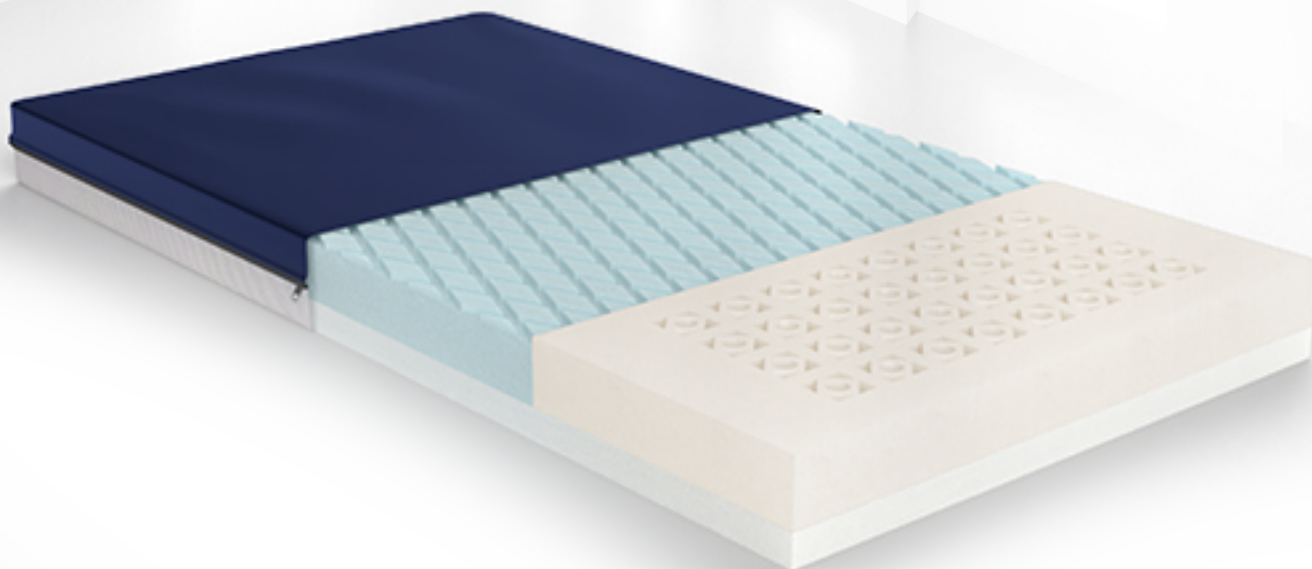

Committed
TO QUALITY

High density, bariatric, therapeutic foam mattress provides pressure redistribution and maximum support.

FEATURES

- Foam heel section for enhanced lower extremity pressure offloading
- Removable premium, multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Durable foam is latex-free and flammability rated CFR 1632/1633

WHAT'S NEW AND IMPROVED?




3 Year Extended
Cover Warranty*



Weight Capacity
Verified



Made In
North America

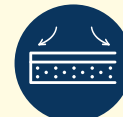

Full Width
Top Layer

Specifications

Model #	1500SC-42-FB	Size	42" (W) X 80" (L) X 6" (H)
Model #	1500SC-48-FB	Size	48" (W) X 80" (L) X 6" (H)
Weight Capacity	650 lbs. (42") - 750 lbs. (48")		
Warranty	5 years limited foam, 3 year cover		


* 20,000 cycle static fatigue testing verified

SUPPORT MATTERS!

NEW!

Gravity 6

LONG TERM CARE THERAPEUTIC FOAM MATTRESS

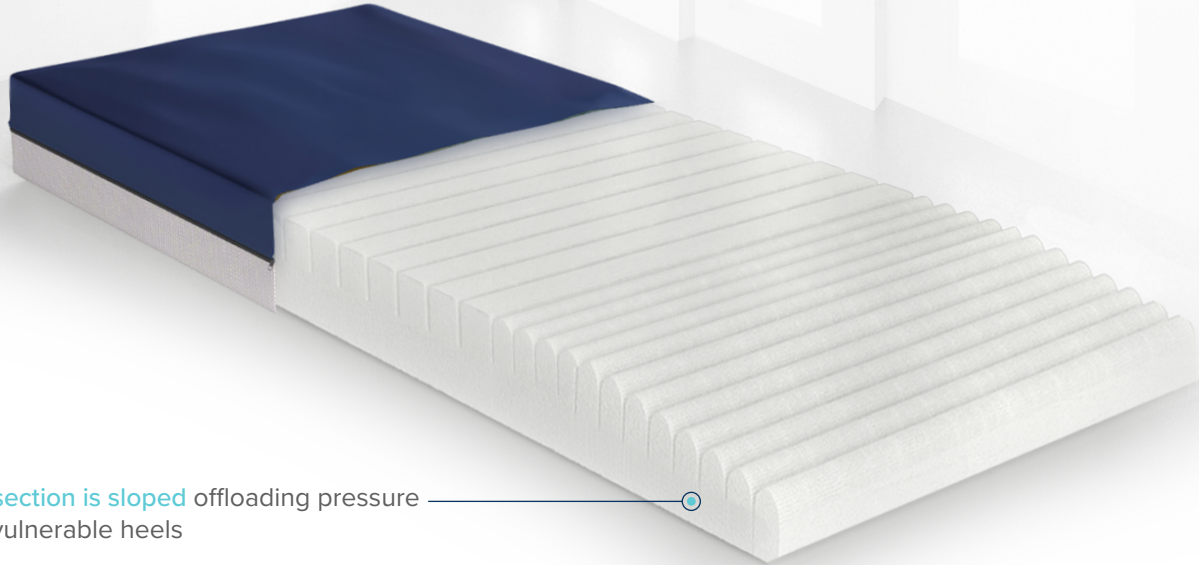

Heel section is sloped offloading pressure from vulnerable heels

ADDITIONAL FEATURES

- Deluxe horizontal, cross-cut foam mattress provides comfort, support, and pressure redistribution over 5 therapeutic pressure zones
- Better compliance and durability than innerspring mattresses allows for improved patient satisfaction and fewer service calls
- Removable multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture.
- Barrier covered zipper minimizes risk of foam contamination
- Non-skid bottom keeps mattress safely in place
- Meets flammability standard CFR 1632 and CFR 1633

Committed
TO QUALITY

High quality latex free therapeutic foam mattress provides pressure redistribution and shear/friction reduction



Barrier Covered
Zipper




3 Year Extended
Cover Warranty*



Weight Capacity
3rd Party Verified



Made In
North America

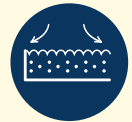

Full Width
Therapeutic Surface

Model Numbers

Model #	15670
Size	36" W x 80" L x 6" H
Weight Capacity	350 lbs.
Warranty	3 years limited foam, 3 years cover

* 20,000 cycle static fatigue testing verified

GROUP 1

Gravity 9

LONG TERM CARE THERAPEUTIC FOAM MATTRESS

NOW AVAILABLE IN 39" WIDTH!
Model #G9-39x80x6

Removable multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture

Full coverage top layer of open cell, visco-elastic memory foam improves immersion and comfort


High resilience foam comfort layer includes high-density, die-cut foam to optimize zoned pressure redistribution

Durable support layer maintains core integrity of the mattress for long lasting individual use

Full length horizontally scored **articulation cut bottom layer** extends the durability and life of the mattress

Channel cut designed **heel slope** offloads pressure to help prevent and treat lower extremity wounds

High density foam perimeter provides solid edge support for patient care and transferring

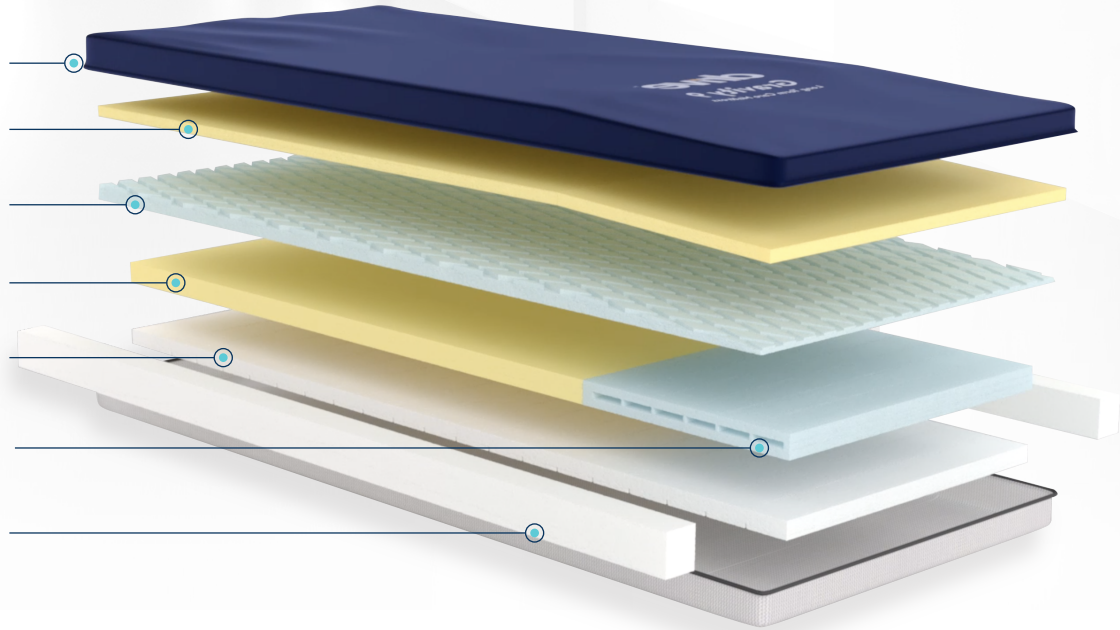

Committed
TO QUALITY

Specifically engineered for pressure redistribution, shear reduction and skin protection, the Gravity series variable density sleep systems include unique combinations of latex-free foam.

ADDITIONAL FEATURES

- Open cell, visco-elastic memory foam full width top layer for pressure redistribution
- 15 year limited foam warranty
- NEW 3 year cover warranty
- 450 lb weight capacity
- CFR 1632/1633 flammability compliant

WHAT'S NEW AND IMPROVED?



3 Year Extended Cover Warranty*




Weight Capacity 3rd Party Verified



New Size Available in 39"



Made In North America

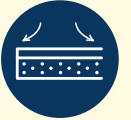

Full Width Top Layer

Model Number	Size	Model Number	Size
15996	36" (W) X 76" (L) X 6" (H)	G9-39x80x6	39" (W) X 80" (L) X 6" (H)
15970	36" (W) X 80" (L) X 6" (H)	15970-4280	42" (W) X 80" (L) X 6" (H)
15984	36" (W) X 84" (L) X 6" (H)	15970-4284	42" (W) X 84" (L) X 6" (H)
15977 (with raised rails)	36" (W) X 80" (L) X 6" (H)		

* 20,000 cycle static fatigue testing verified